

# Particle View Software

## Particle Counter Download Utility

---

Particle View software includes all of the functions required for the downloading, viewing of data, and production of both tabular and graphic reports. It is designed for use with the Met One Instruments, Inc. line of hand held and bench top particle counters and is compatible with the following models: GT-324, GT-521, GT-521S, GT-526, GT-526S, 804, BT-637, BT-637S, BT-610 and BT-620.

### Particle View Features:

#### Passwords

Two password levels (administrator and user).

#### Secure Database

Windows SQL Server Express database.

#### Unit Conversions

Concentration conversions (particles / liter, particles / cubic foot, particles / cubic meter) and temperature (°C or °F).

#### Alarm Limits

Particle count, temperature and relative humidity alarm limits (High, High-High, Low, Low-Low).

#### Sort / Filter Data

Sort data by date / time, location and particle counts. Filter data by date / time, location, and particle size.

#### Reports

View and print tabular reports. Export spreadsheet compatible files (CSV format).

#### Charts

View, print and export graphical charts.

### System Requirements:

Processor Speed:	1GHz or faster
Memory:	512MB minimum
Disk Space:	60MB
Operating System:	Windows XP, Vista, Windows 7, Windows 8 and Windows 10

### Applications:

Cleanrooms	Indoor Air Quality	Filter Testing	Emission Sourcing
Controlled Environments	Process Control	Compliance Testing	Remote Sampling




# Particle View Images

## Main Display

Date/Time	Loc ID	0.3u	0.5u	0.7u	1.0u	2.0u	5.0u
4/23/2012 04:33:34	001	139,770	16,060	7,050	4,880	2,170	340
4/23/2012 04:32:34	001	141,120	15,840	7,160	5,020	2,290	530
4/23/2012 04:31:34	001	144,300	16,890	8,020	5,460	2,510	630
4/23/2012 04:30:34	001	148,540	17,260	8,140	5,570	2,770	850
4/23/2012 04:29:34	001	148,440	18,360	8,830	6,550	3,380	1,120
4/23/2012 04:24:14	001	154,250	14,750	6,630	4,570	2,310	640
4/23/2012 04:23:14	001	159,180	15,360	6,690	4,650	2,100	500
4/23/2012 04:22:14	001	150,630	15,290	6,740	4,800	2,580	1,110
4/23/2012 04:21:14	001	140,030	12,930	5,540	3,680	1,810	600
4/23/2012 04:20:14	001	139,010	13,010	5,070	3,250	1,630	830

## Chart Display

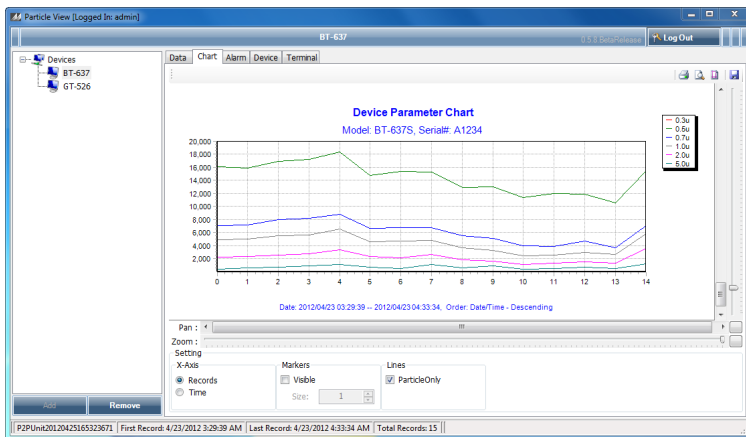

## Report Preview

Date/Time	Loc	0.3u	0.5u	0.7u	1.0u	2.0u	5.0u
2012/04/23 04:33:34	001	139,770	16,060	7,050	4,880	2,170	340
2012/04/23 04:32:34	001	141,120	15,840	7,160	5,020	2,290	530
2012/04/23 04:31:34	001	144,300	16,890	8,020	5,460	2,510	630
2012/04/23 04:30:34	001	148,540	17,260	8,140	5,570	2,770	850
2012/04/23 04:29:34	001	148,440	18,360	8,830	6,550	3,380	1,120
2012/04/23 04:24:14	001	154,250	14,750	6,630	4,570	2,310	640
2012/04/23 04:23:14	001	159,180	15,360	6,690	4,650	2,100	500
2012/04/23 04:22:14	001	150,630	15,290	6,740	4,800	2,580	1,110
2012/04/23 04:21:14	001	140,030	12,930	5,540	3,680	1,810	600
2012/04/23 04:20:14	001	139,010	13,010	5,070	3,250	1,630	830

Specifications subject to change without notice.

Rev. January 2020



1600 Washington Blvd. • Grants Pass, OR 97526 • 541.471.7111 • www.metone.com  
[sales@metone.com](mailto:sales@metone.com)